



Designing the Next Generation Liquid Nitrogen Biological Repository

Innovation. Experience. Performance.®



Traditional Temperature Modalities

- Ambient
 - Temperature / Humidity Control
- +4°C, -20°C
 - Mechanical refrigeration
 - ❖ Walk-in
 - ❖ Reach-in refrigerator
 - Very reliable, least complex mechanical refrigeration system
- -80°C
 - Mechanical refrigeration
 - ❖ Walk-in
 - ❖ Chest or upright freezer
 - Increasingly complex, decreasing reliability
- < -135°C
 - Mechanical refrigeration, highly complex and prone to mechanical failure
 - Liquid Nitrogen Refrigeration, least complex, most reliable refrigeration system



Common Components of Biorepository

- Secure facility
- Emergency Backup power
- Enhanced HVAC capacity
- Centralized alarm system
- Bulk LN2 and piping system
- O2 monitoring
- Freezers...



Infrastructure to Support a Modern Biorepository

ULT Storage


■ ULT Storage

- Complex electrical distribution system
 - ❖ Emergency power generation
 - ❖ Line conditioning
 - ❖ High voltage / high amperage capacity to building
- Highly enhanced HVAC system
 - ❖ Multiple HVAC systems added based on number of mechanical freezers
 - ❖ Complicated duct work to evenly distribute cooling capacity
- Significant backup specimen storage
 - ❖ Capacity requirement increases with average freezer age
- Highly qualified refrigeration service technician, in-house or contracted



Infrastructure to Support a Modern Biorepository

ULT Storage

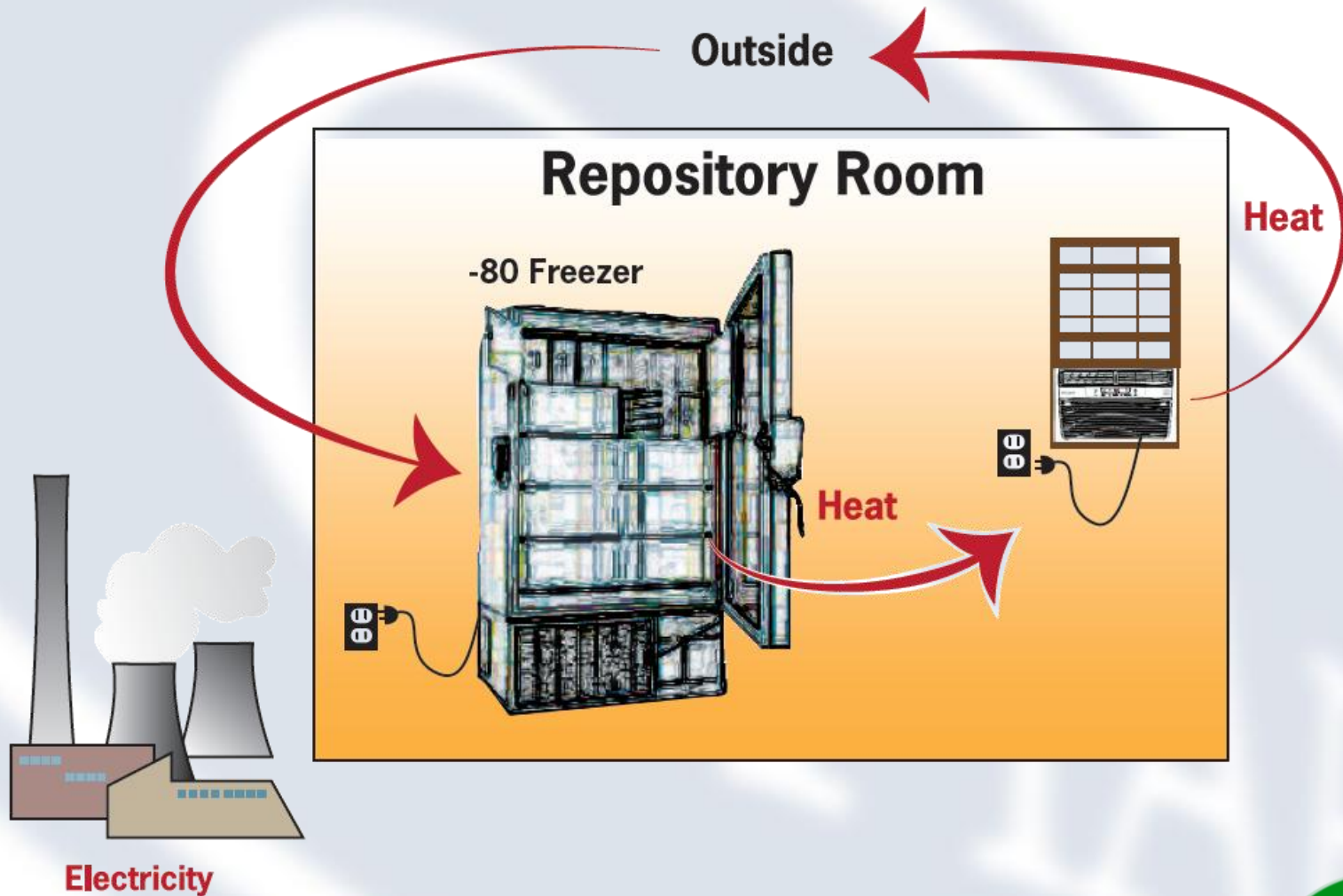

Infrastructure to Support a Modern Biorepository


Cryogenic Storage

- Liquid nitrogen source, Dewars or bulk distribution
- Oxygen depletion alarms
- Standard electrical installation with optional backup power
 - Household power requirements
 - Units available with battery backup to limit backup power requirement
- Standard HVAC system with enhanced ventilation capacity
 - Most commercial HVAC systems adequate for installation



Infrastructure to Support a Modern Biorepository

Cryogenic Storage

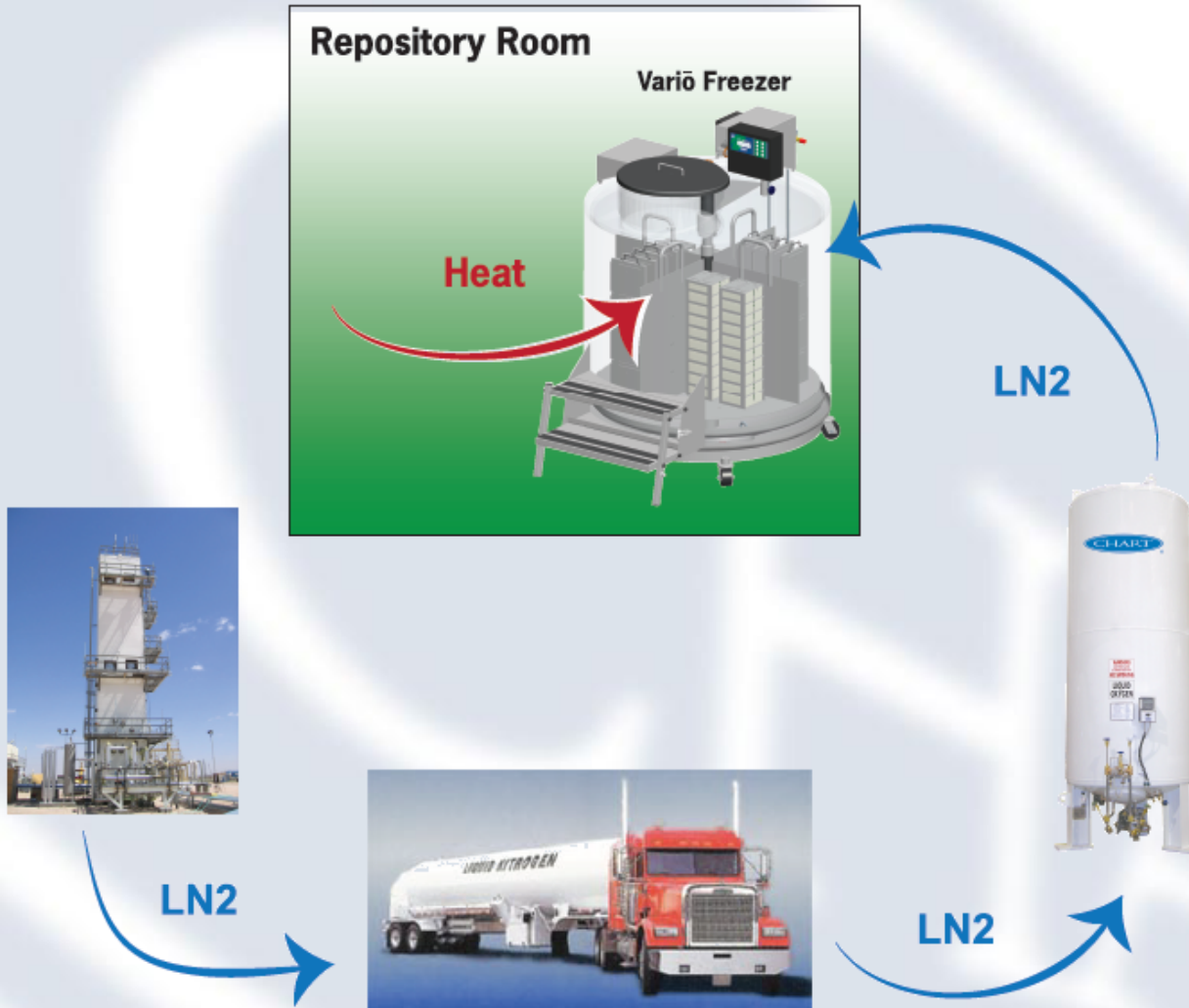


ULT vs Cryogenic

Complexity Scale



Increasing Complexity

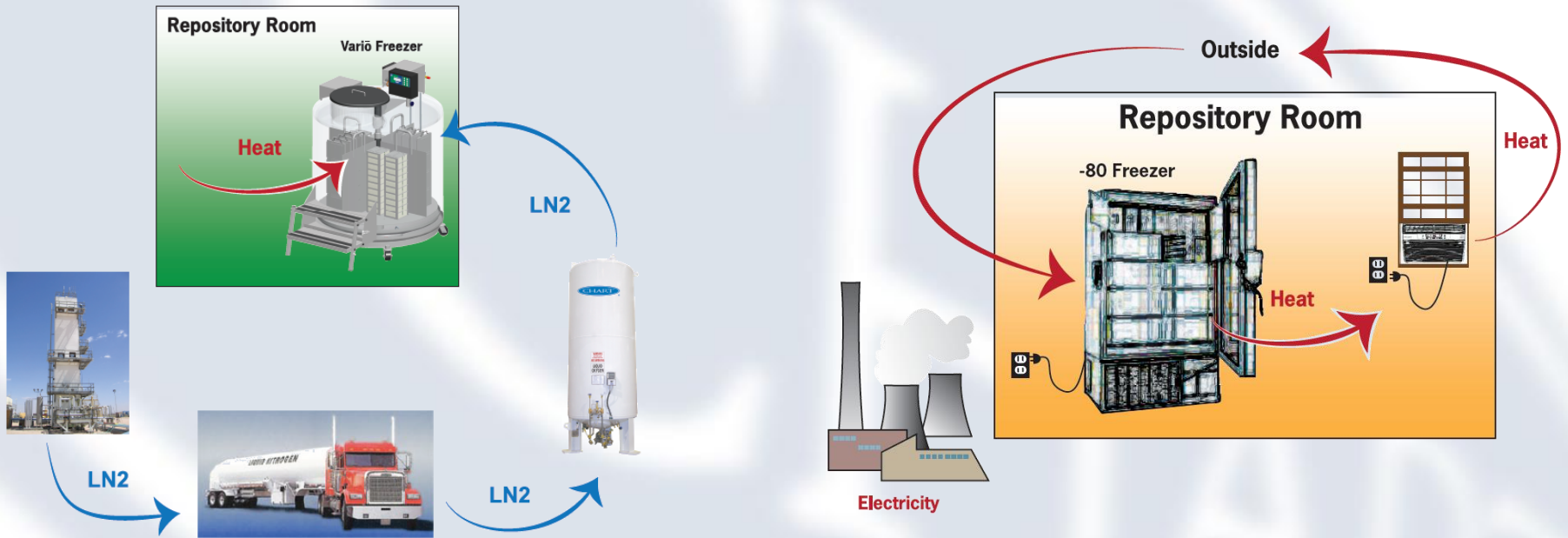

What Keeps You Up at Night?

- Power loss
- False alarms
- Everyday something breaks...
 - Compressor failure
 - Electronic failure
 - HVAC failure
 - You name it: _____

90% of your headaches arise from mechanical refrigeration



Reasons LN2 storage helps you sleep easy

- Simplicity of LN2 has made storage of material at cryogenic temperatures a highly reliable segment of a biorepository facility
 - Long freezer hold time increases your safety margin and response time window
 - Bulk storage of LN2 guarantees on-site refrigeration
 - ❖ LN2 supply highly reliable
 - Minimal moving parts assures a highly reliable and robust system
 - ❖ No compressors: freezers or HVAC system



Chart MVE and Cryo Associates Introduce the Reliability of LN2 at ULT

- For years we tried to convince you that colder is better
- With great pain, we have accepted that the -80°C freezer is here to stay.
- With great excitement, we are introducing the Vario Freezer.





Variō Series



Innovation. Experience. Performance.®



MVE Variō Design

- Innovative and energy efficient alternative for ultra-low temperature to cryogenic storage
- User-defined temperature anywhere between -50°C and -150°C
- Dry sample space
 - Significantly reduces risk of sample contamination via contact with LN₂
- Consistent temperature profile
 - Even with the lid open
- Extended hold time safety margin
- Less than 1% of power consumption compared to leading upright mechanical freezers

Low maintenance and long lasting
doubled-walled, vacuum insulated,
all-stainless steel construction




Advanced MVE Variō™ Pro controller

Interior heat exchange system: LN₂
flows through heat exchangers and
the vaporization energy provides
low-cost refrigeration





Innovative, patent-pending refrigeration system

- LN₂ flows through a heat exchange system inside the freezer
- Vaporization energy of LN₂ cools the freezer
- Fully utilizes the heat capacity of LN₂ while purging moisture from the system
- LN₂ metered and monitored to maintain user-defined temperature

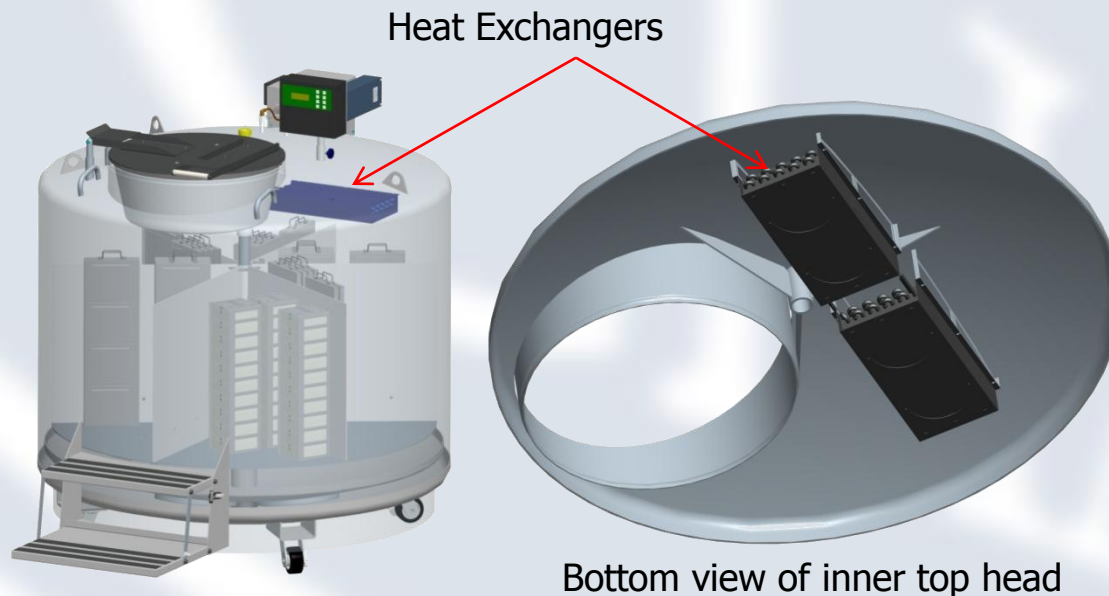

Cooling Cycle


LN₂ begins to flow into the heat exchangers inside the freezer. The vaporization energy of the LN₂ cools the freezer until the low end of the temperature range is reached.

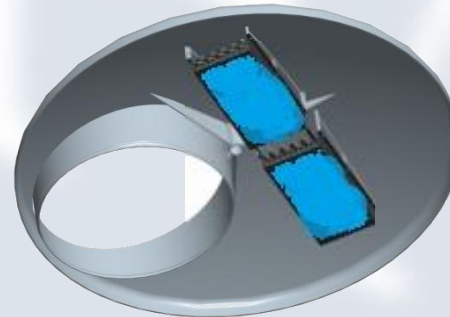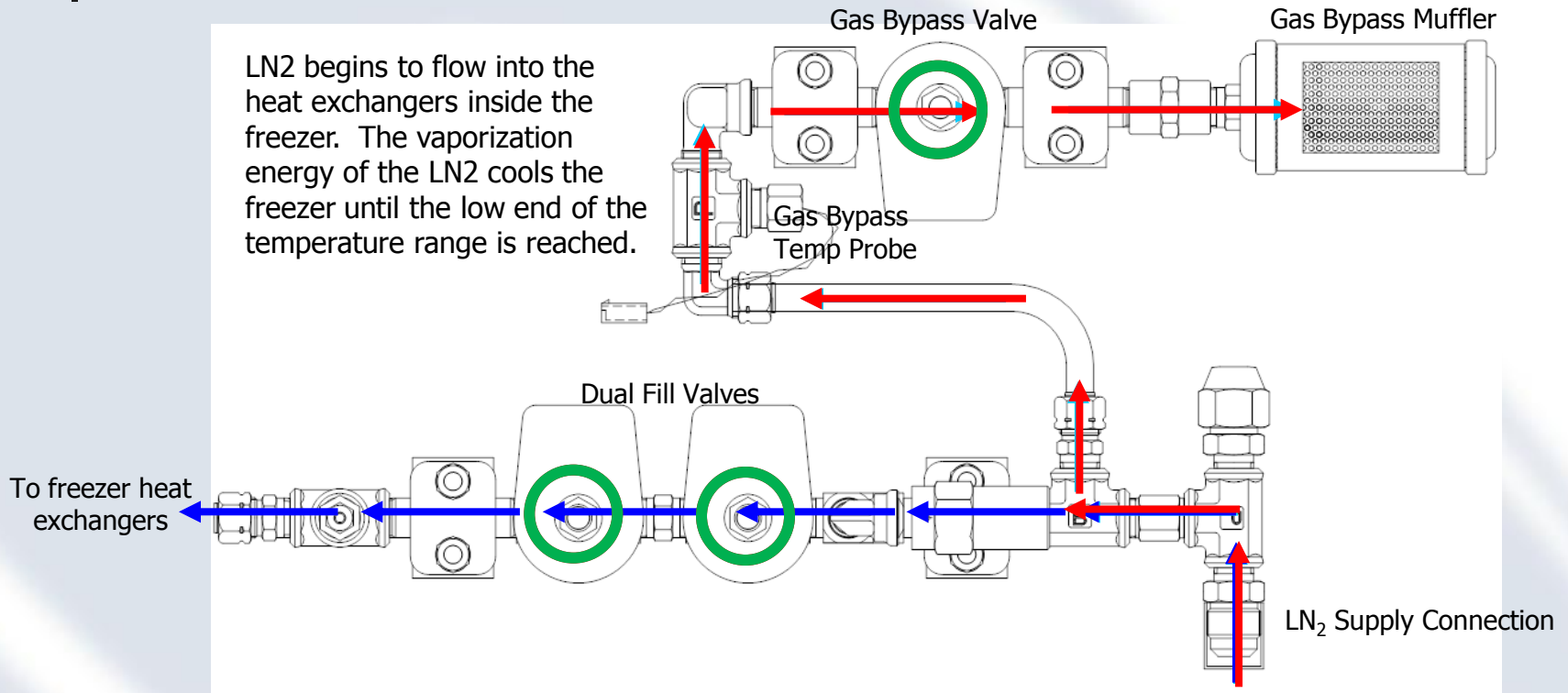

Temperature Profile -80°C


Consistent temperature profile throughout the storage space, even with the lid open.

While the lid is open, the cooling cycle interval is compressed, but the storage space temperature does not increase outside of the set range.

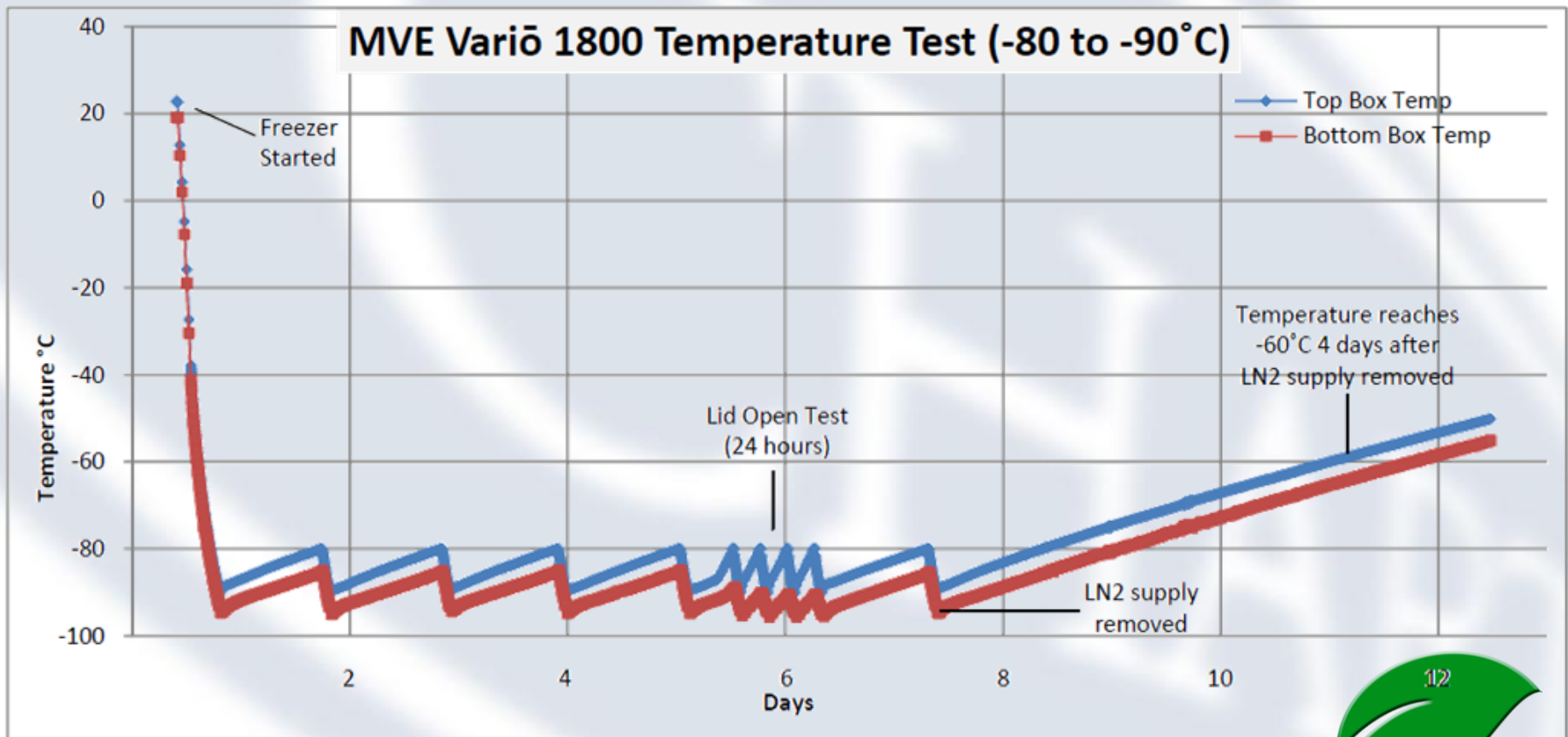


Temperature Gradient

Temperature gradient within storage space (top box to bottom box) is between 2°C and 5°C

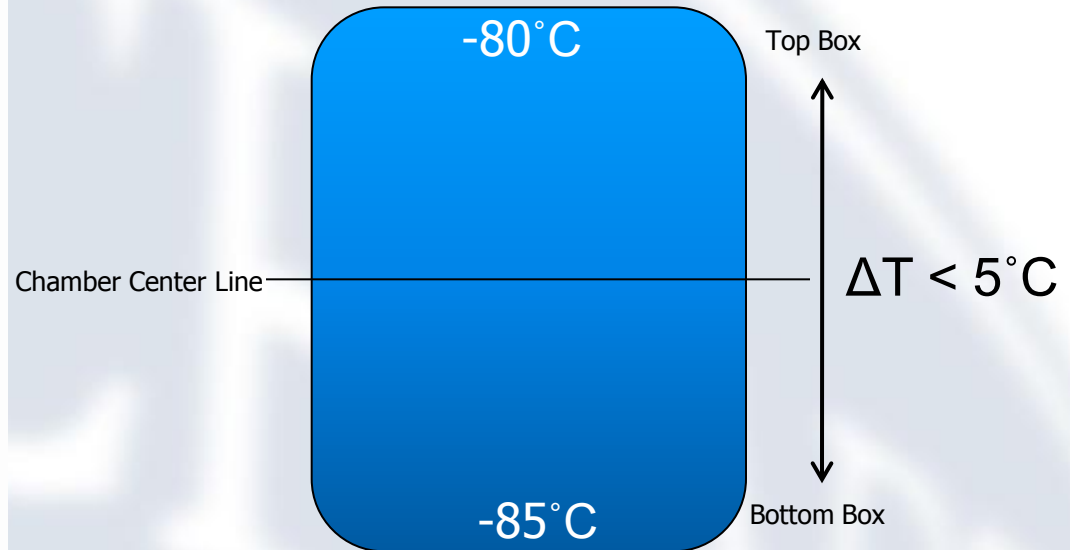

Lid Open Security

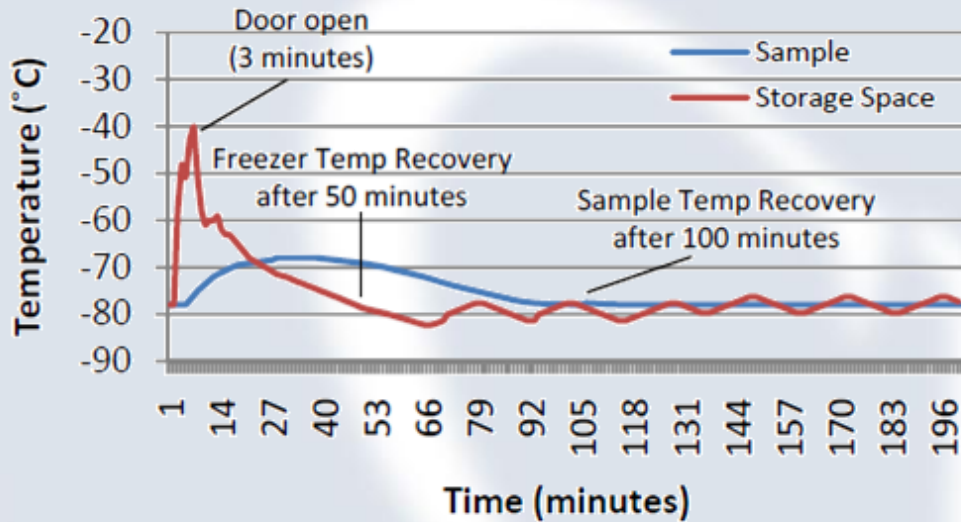

- Typical upright mechanical freezer
- Significant temp change and recovery time


Does not compare to MVE's sample security



Hold Time Safety Margin

- Refrigeration system failure until critical temperature is reached
- The MVE Variō™ 1800 Series provides by far the longest hold time in the industry.
- When the LN2 supply is depleted, it takes 4 days (96 hours) when storing at -80°C to rise above -60°C


Hold Time Comparison from -80°C to -60°C

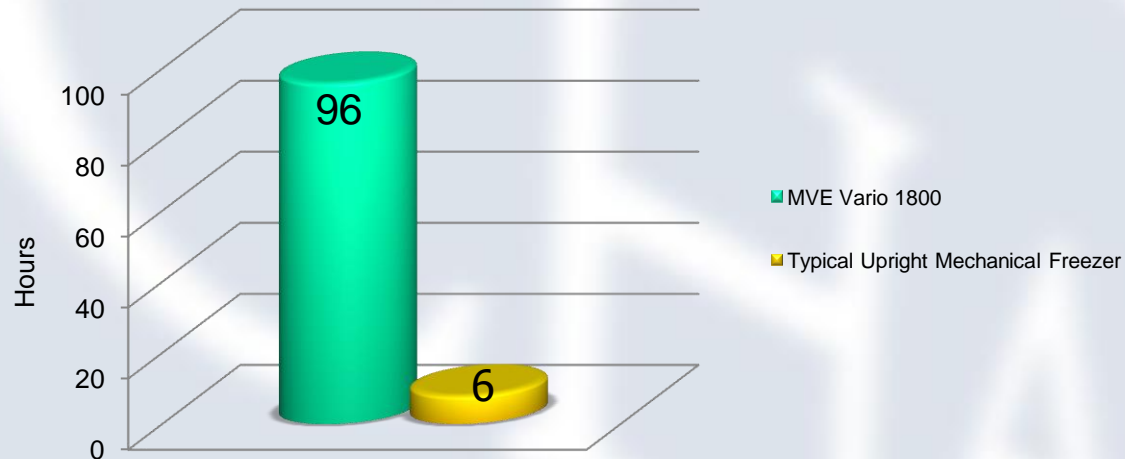

Power Consumption


Up to 81,900 2mL vials while consuming less than 1% of the power compared to the leading upright mechanical freezers, only 8 W continuous.

Does not account for the additional HVAC requirements for mechanical freezers.

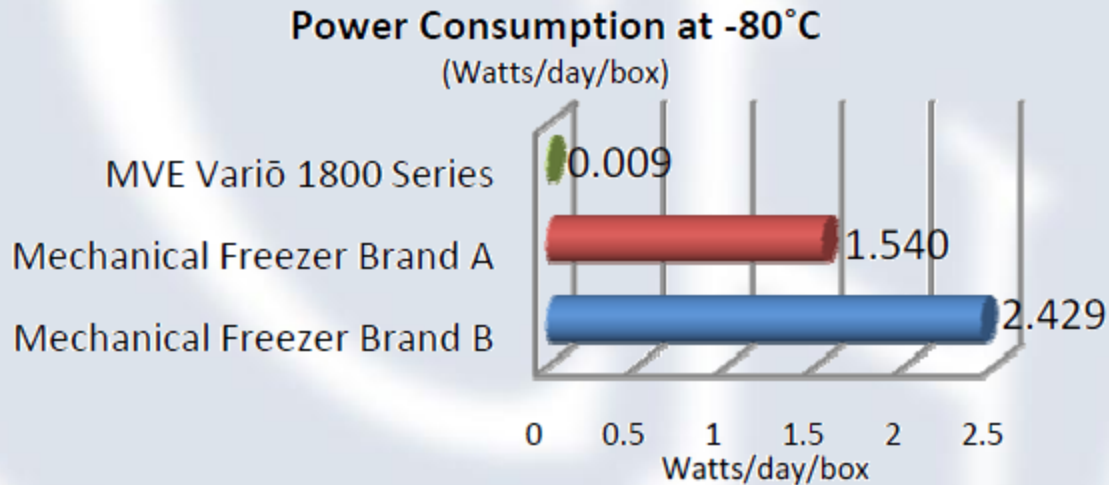

Insulation Deterioration over time

■ MVE Vacuum Insulation System ■ Traditional Mechanical Freezer

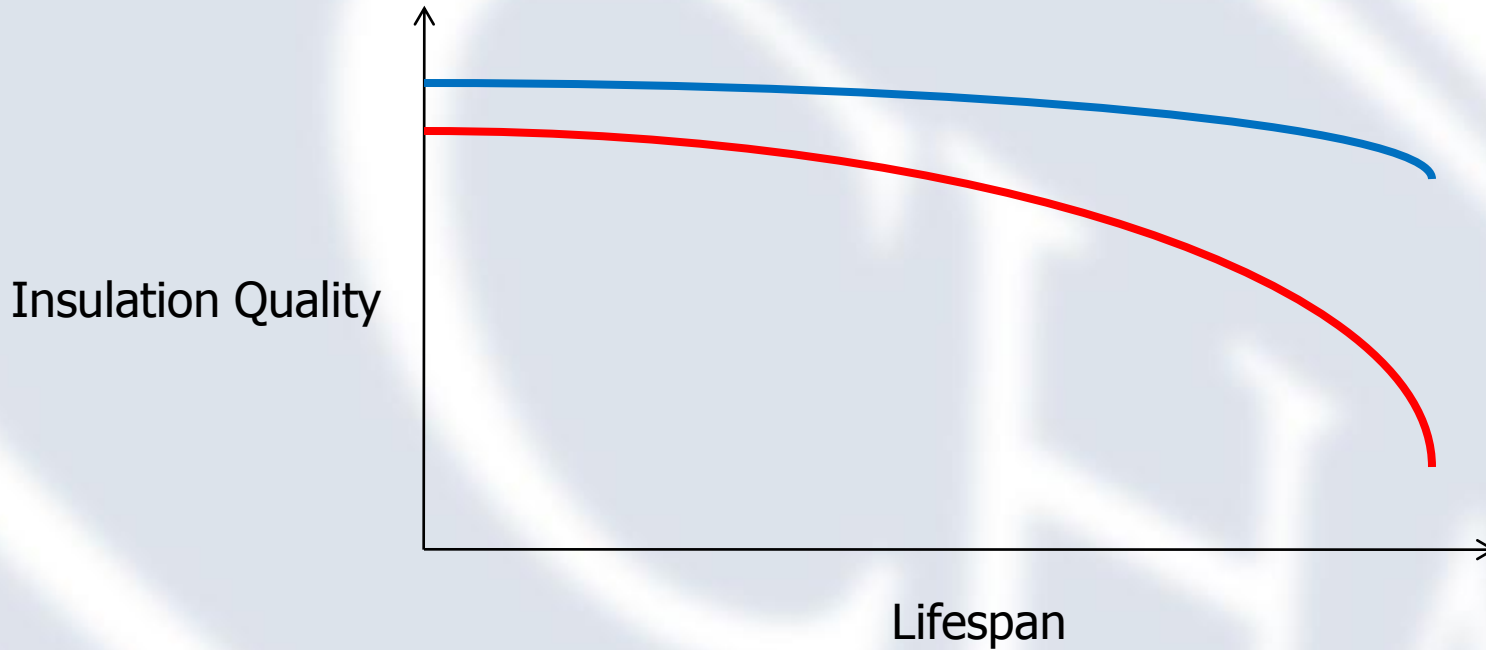



Summary

- Highest sample security
 - Lid open temperature security
 - Hold Time Safety Margin
 - Reduced risk of LN2 contamination
- ROI Potential
 - Less than 1% of the power consumption compared to mechanical freezers
 - Significantly reduced operating costs
 - Highest sample storage density
- Energy efficient and environmentally friendly alternative
 - Minimal power consumption – similar to that of an alarm clock
 - No additional HVAC requirements
 - Zero CO2 emissions
 - No ozone depleting refrigerants; CFC or HFC
 - No disposal issues – over 90% recyclable



Case Study: 4M Sample Repository

- Chart and Cryo Associates surveyed several repositories over the past two years to understand the Capital Investment and Operating Costs of a repository.
- From the data collected, we created cost models for two different types of repositories, both with a 4M sample capacity
 - Repository A:
 - Mechanical -80 Capacity: 75%
 - Liquid Nitrogen Capacity: 25%
 - Repository B:
 - Vario -80 Capacity: 75%
 - Liquid Nitrogen Capacity: 25%
- Vario -80 based repository shows a substantial reduction in operating costs over traditional -80 mechanical storage



Repository A: 75% Upright Mechanical

	Primary	Backup	Total Freezers	Samples/Freezer	Total Samples
LN2 Freezers	20	2	22	67,080	1,341,600
Upright -80	70	7	77	38,880	2,721,600
				Grand Total:	4,063,200

Assumptions:

- 81 cell boxes
- Standard size rack use only – no “Mini Racks”
- Total sample capacity does not include backup freezers



Repository A: Capitol Costs

Capital Costs	LN2	Upright -80	Initial Investment	Additional Cost Year 10	Additional Cost Year 20
Freezer Cost (including racks)	\$38,000	\$9,962	\$1,603,074	\$280,544	\$359,117
Bulk Tank / VIP 3000 gallon	\$178,000		\$178,000	\$0	\$0
Alarm System	\$2,000	\$1,000	\$121,000	\$0	\$121,000
Backup Power (250 kW)		\$80,000	\$80,000	\$0	\$80,000
HVAC			\$176,000	\$0	\$140,000
O2 Monitoring			\$6,000	\$0	\$6,000
			\$2,164,074	\$280,544	\$706,117
			Total Capital Investment over 30 years		\$3,150,735

Assumptions:

- Cost analysis based on freezer storage components
- Freezer costs include racking systems
- Upright Mechanical Replacement every 10 years
- Annual price increase for replacements 2.5%
- 116 Ton HVAC system priced at \$1,500/ton installed, \$1200/ton replacement costs



Repository A: Operating Costs

Operating Costs	LN2	Upright -80	1 Year	10 Years	20 Years	30 Years
Electricity (including HVAC)	1	75	\$69,564	\$695,640	\$1,391,280	\$2,086,920
LN2	60		\$15,840	\$158,400	\$316,800	\$475,200
Generator PM	0	2,800	\$2,500	\$25,000	\$50,000	\$75,000
Freezer PM	135	85	\$9,515	\$95,150	\$190,300	\$285,450
Freezer Repair	190	380	\$33,440	\$334,400	\$668,800	\$1,003,200
Freezer Labor	115	725	\$58,355	\$583,550	\$1,167,100	\$1,750,650
			\$189,214	\$1,892,140	\$3,784,280	\$5,676,420
				Total Operating + Capital		\$8,827,155

Assumptions:

- Includes re-vac of freezer every 10 years
- Includes periodic de-icing/defrosting of freezers
- Includes daily performance check of all freezers
- Does not include after hour/emergency repairs



Repository A: Summary

	At Start-up	After 30 Years
Capitol Costs	\$2,100,000	\$3,100,000
Operating Costs	\$190,000	\$5,700,000
Total:	\$2,190,000	\$8,800,000



Repository B: 75% Vario N₂ Vapor

	Primary	Backup	Total Freezers	Samples/Freezer	Total Samples
LN2 Freezers	20	2	22	67,080	1,341,600
Vario 1881R	40	4	44	67,080	2,683,200
				Grand Total:	4,024,800

Assumptions:

- 81 cell boxes
- Standard size rack use only – no “Mini Racks”
- Total sample capacity does not include backup freezers



Repository B: Capitol Costs

Capital Costs	LN2	Vario 1881R	Intial Investment	Additional Cost Year 10	Additional Cost Year 20
Freezer Cost (including racks)	\$38,000	\$40,000	\$2,596,000	\$0	\$0
Bulk Tank / VIP (9000 gallon)			\$320,000	\$0	\$0
Alarm System	\$2,000	\$1,000	\$88,000	\$0	\$88,000
Backup Power (250 kW)	\$0	\$0	\$8,000	\$0	\$8,000
HVAC	\$0	\$0	\$17,600	\$0	\$14,000
O2 Monitoring			\$6,000	\$0	\$6,000
			\$3,035,600	\$0	\$116,000
			Total Capital Investment over 30 years		\$3,151,600

Assumptions:

- Freezer costs include racking systems
- HVAC costs approximately 10% of requirement for Mechanical



Repository B: Operating Costs

Operating Costs	LN2	Vario 1881R	1 Year	10 Years	20 Years	30 Years
Electricity (including HVAC)	\$1	\$1	\$792	\$7,920	\$15,840	\$23,760
LN2	\$60	\$60	\$47,520	\$475,200	\$950,400	\$1,425,600
Generator PM	\$0	\$0	\$0	\$0	\$0	\$0
Freezer PM	\$135	\$135	\$8,910	\$89,100	\$178,200	\$267,300
Freezer Repair	\$190	\$190	\$12,540	\$125,400	\$250,800	\$376,200
Freezer Labor	\$115	\$115	\$7,590	\$75,900	\$151,800	\$227,700
			\$77,352	\$773,520	\$1,547,040	\$2,320,560
					Total Operating + Capital	\$5,472,160

Assumptions:

- Includes re-vac of freezer every 10 years
- Includes daily performance check of all freezers
- Does not include after hour/emergency repairs



Repository B: Summary

	At Start-up	After 30 Years
Capitol Costs	\$3,100,000	\$3,100,000
Operating Costs	\$77,000	\$2,300,000
Total:	\$3,177,000	\$5,400,000



Comparison Summary

	Repository A	Repository B	Difference
Initial Capital Costs	\$2,100,000	\$3,000,000	+\$900,000
Capital Costs after 30 years	\$3,100,000	\$3,100,000	-
1 st year operating costs	\$189,000	\$77,000	-\$112,000
Operating Costs 30 years	\$5,700,000	\$2,300,000	-\$3,400,000
Operating + Capital 30 yrs	\$8,800,000	\$5,400,000	-\$3,400,000

Assumptions

- Price of peace of mind not factored into analysis

